


ATAC-seq Data Pre-processing Pipeline

Chao-Jen Wong
October 1, 2018

Overview

Buenrostro 2015:

Assay for Transposase Accessible Chromatin with high-throughput sequencing (ATAC-seq) is a method for mapping chromatin accessibility genome-wide.

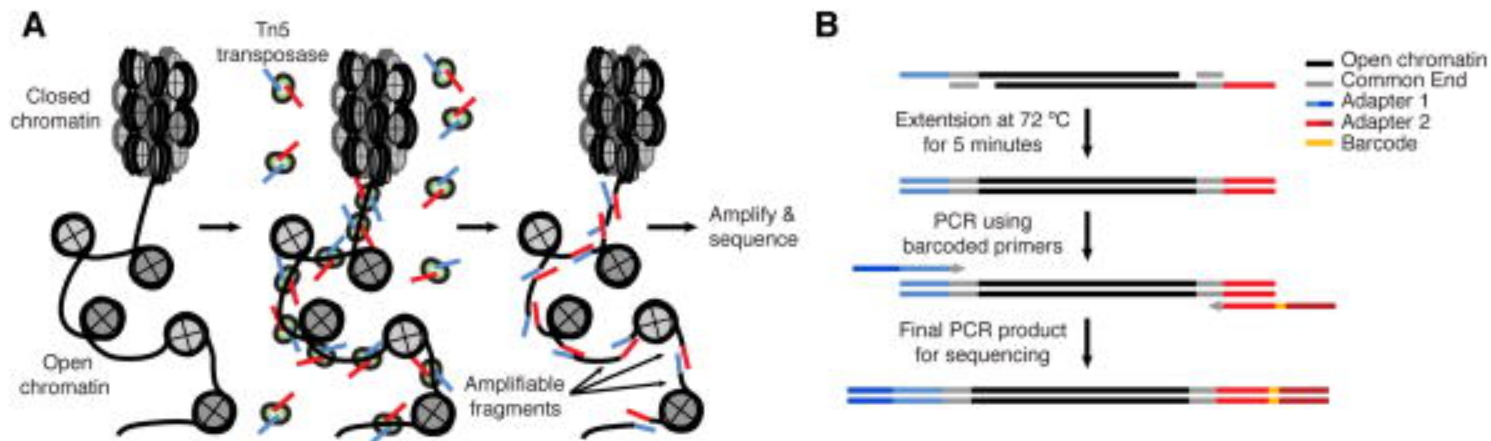

Data analysis

Data Pre-processing: make clean BAM files

- Adaptor removal
- Alignments
- PCR duplicates
- Exclude chrM
- Sanity check

Downstream analysis

- Peaking calling
- Annotation
- Differential binding

Data analysis

Software requirement

- cutadapt
- picard
- samtools
- bowtie2
- FASTQC
- MACS2

BioConductor

- ChIPseeker or ChIPpeakAnno
- DiffBind

Pre Pre-Processing

Quality Control

FASTQC

Chasity filtering

Casava 1.8 format -> retain reads with filter flag=N
@EAS139:136:FC706VJ:2:2104:15343:197393 1:N:18:ATCACG

```
for i in *fastq.gz
do
    i2=${i//.gz/}
    zgrep -A 3 '^@.*[^:]*:N:[^:]*:' $i \
    | zgrep -v '^\\-\\-$' > $fileredDir/$i2
done
```

Tips: Huthbase.org offers FASTQC+Filtering pipeline

Pre Pre-Processing

Concatenate sequence reads on multiple lanes

Genomic Core sequence samples on multiple lanes may have separate files for each lane.

```
cat $(ls *_R1_*.fastq.gz) > sampleName_R1.fastq.gz
```

or

```
zcat $(ls *_R1_*.fastq.gz) > sampleName_R1.fastq
```

Adapter removal

[cutadapt](#) searches input reads for a given adapter sequence and removes it.

Usual parameters

- `-a` to remove 3' adapter
- `-g` to remove 5' adapter
- `-o` to give minimal length of the provided adapter sequence.

Adapter removal

Illumina Nextera Transposase Adapters

Read 1

5' TCGTCGGCAGCGTCAGATGTGTATAAGAGACAG

Read 2

5' GTCTCGTGGGCTCGGAGATGTGTATAAGAGACAG

```
ml cutadapt
```

```
cutadapt -g TCGTCGGCAGCGTCAGATGTGTATAAGAGACAG \  
  -O 5 -e 0.2 \  
  input_R1.fastq > output_R1.fastq 2> report_R1.txt  
gzip output.fastq
```


Alignment using Bowtie2

Genomic indexes

Pre-build reference are available from iGenome

```
cd /shared/biodata/ngs/Reference/iGenomes
cwon2@rhino2:/shared/biodata/ngs/Reference/iGenomes$ ls -lt
total 682
drwxrwsr-x 3 rbasom biodata_adm 22 Dec 6 2017 Bos_taurus
drwxrwsr-x 7 mfitzgib biodata_adm 155 Jul 27 2017 Homo_sapiens
drwxrwsr-x 4 rbasom biodata_adm 47 Aug 22 2016 Rattus_norvegicus
drwxrwsr-x 3 rbasom biodata_adm 25 Jun 8 2016 Macaca_mulatta
drwxrwsr-x 4 rbasom biodata_adm 47 Jan 26 2016 Caenorhabditis_elegans
drwxrwsr-x 4 rbasom biodata_adm 47 Jan 14 2016 Saccharomyces_cerevisiae
drwxrwsr-x 3 rbasom biodata_adm 22 Nov 4 2015 Danio_rerio
drwxrwsr-x 3 rbasom biodata_adm 22 Sep 1 2015 Gallus_gallus
drwxrwsr-x 4 rbasom biodata_adm 47 May 18 2015 Canis_familiaris
drwxrwsr-x 4 mfitzgib biodata_adm 44 Oct 30 2014 Drosophila_melanogaster
```

Alignment using Bowtie2

Usual paramters

parameter	description	default
-x <int>	Maximum DNA fragment length	500
--very-sensitive	better alignment results are frequently achieved with it	-sensitive
-p <int>	Number of cores on which to run	1
--seed <int>	number used to initialize a pseudorandom number generator	
-k <int>	Search for one or more alignments, report alignments	
-a	Search for and report all alignments	

Alignment using Bowtie2

```
ml bowtie2/2.2.5 # if not already loaded
ml samtools/1.0
bowtie2 --very-sensitive -p 4 \
    -I 20 -X 1000 -x $hg38_genomeBuild \
    -1 $trim1.gz -2 $trim2.gz \
    | samtools view -Sb - > $sampleName.bam
```

Alignment using Bowtie2

Why setting **-x 1000** is reasonable?

- nucleosome (~147bps) + space in between (~20bps) ~= 160 bps
- cover up at least three size of nucleosomes ~= 600 bps
- cover up possible skipped nucleosomes near TSS
- just to loosen up

Sort BAM

```
# sort  
samtools sort -@ 4 $sampleName.bam $sampleName.bam.sorted  
mv $sampleName.bam.sorted.bam $sampleName.bam
```

Exclude chrM

ATAC-seq datasets usually contain decent amount of reads derived from chrM.

```
samtools view -h $sampleName.bam | awk '{if($3 != "chrM"){print $0}}' \
| samtools view -Sb - > $sampleName\_rmChrm.bam
```

```
mv $sampleName\_rmChrm.bam $sampleName.bam
```

Very slow. Better way: write a python package or use picard(?) to remove chrM.

PCR duplicates

PCR duplicates are exact copy of fragments and are artifacts of the library preparation.

```
ml picard/2.7.1-Java-1.8.0_92
java -jar ${EBROOTPICARD}/build/libs/picard.jar MarkDuplicates \
    I=$sampleName.bam O=$sampleName\_filter.bam \
    M=$sampleName\_dup\_metrics.txt REMOVE_DUPLICATES=true \
    ASSUME_SORTED=true \
    TMP_DIR=$HOME/tmp

mv $sampleName\_filer.bam $sampleName.bam
```

Tips: use sorted BAM and TMP_DIR to save space and time

Index BAM

Finish up!

```
samtools index $sampleName.bam  
touch $sampleName.bowtie2Done.txt
```


Peak Calling using MACS2

Example of using proper paired reads and real insert size:

```
ml MACS2/2.1.1.20160309-foss-2015b-Python-2.7.9
macs2 callpeak -f BAMPE -t $sampleName.bam \
  -g hs \
  --outdir <output_dir> -n <file_prefix> -q 0.05 \
  --keep-dup all
```

Peak Calling using MACS2

Hints

- If PCR duplicates are not removed already, use default setting `--keep-dup 1` to keep one tag of at the sample position.
- If looking for enriched cutting sites. Assuming 200 is the average fragment size:

```
--nomodel --shift -100 --extsize 200
```

Annotation

Bioconductor **ChIPseeker** or **ChIPpeakAnno** packages

```
library(ChIPseeker)
library(plyranges)
library(ggplot2)
library(org.Hs.eg.db)
library(TxDb.Hsapiens.UCSC.hg38.knownGene)
```

Annotaion

Bridge MACS2's narrowPeaks bed file to ChIPseeker

```
#' peaks bed file -> GRanges
readNarrowPeaksBed <- function(bed_file) {
  peaks <- read.table(bed_file, sep="\t", stringsAsFactors=FALSE)
  colnames(peaks) <- c("seqnames", "start", "end",
                      "name", "int10Log10qvalue", ".",
                      "fold_enrichement", "nLog10pvalue",
                      "nLog10qvalue", "summit")
  plyranges::as_granges(peaks)
}
```

Annotation

Example using `ChIPseeker::annotatePeak()` and `ChIPseeker::plotAnnoBar()`.

```
x <- "sampleName_peak.narrowPeaks"
peaks <- readNarrowPeaksBed(x)
peaks_annotation <-
  ChIPseeker::annotatePeak(peaks,
                           TxDb = TxDb.Hsapiens.UCSC.hg38.knownGene,
                           annoDb = "org.Hs.eg.db")
gg <- ChIPseeker::plotAnnoBar(peaks_annotation)
ggsave(filename=filename, gg, device="pdf", width=6, height=3)
```

Annotation

Example using `ChIPseeker::annotatePeak()` and `ChIPseeker::plotAnnoBar()`

